




Other dental-related products

With over 100 years of commitment to quality, we prioritise safety and reliability. To meet the diverse needs of our customers, we offer a comprehensive range of products.

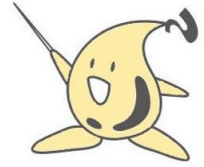

Separator for dental resin bases

Mysep

Mysep is a denture base resin separator with sodium alginate as its main component. It effectively works as a separator between gypsum and resin, and between gypsum and flask.



500g plastic container

Separator for dental resin bases

Features

- Smooth Application: Allows for smooth application, forming a strong film on the gypsum surface.
- Easy Separation: Facilitates easy separation of the resin base, resulting in a clean resin surface.

Gypsum dissolving solution

Quick Melty

Quick Melty can remove gypsum and gypsum-based investment materials attached to dentures and castings by simply soaking them. It dissolves gypsum at approximately twice the speed of our previous products.



1000ml plastic bottle

Dissolving and removing gypsum attached to resin-based dentures
Dissolving and removing gypsum-based investment materials attached to castings

Features

- Simple Just soak in the undiluted solution. When the bubbling stops, it's time to replace the solution.
- Fast The dissolution time for 1g of hardened gypsum is 26 minutes for hard gypsum and 14 minutes for dental calcined gypsum.

Powder storage container

Multi-Powder Box

The Multi-Powder Box features a one-touch open lid. Its sealed design makes it suitable for storing gypsum and impression materials.



1 piece per box

Storage of various gypsum and powder impression materials.

Features

- One-Touch Open The lid opens with a light touch of the button in the center of the lid.
- Easy to Attach and Detach The lid and scraping board frame can be easily attached and detached, making it easy to refill gypsum or impression materials and replace aluminum packs. It can also be washed thoroughly.
- Easy to Retrieve The retrieval opening is tilted forward, making it easier to retrieve materials even from a high position.
- Non-Slip Base The bottom of the box has a non-slip base, providing stability and preventing the box from moving when scooping powder.

For inquiries and product orders, please contact the following



Manufacturer and Distributor

SAN ESU GYPSUM CO.,LTD

<https://www.sanesugypsum.com/dental/>